Therabody

SmartGoggles™

Powered by SmartSense Technology™

User Manual

Table of contents

EN	3-10
ES	11-15
FR	16-20
IT	21-25
DE	26-30
ESMX	31-35
FRCA	36-40
SC	41-45
TC	46-50
KR	51-55
JP	56-60

Born in Los Angeles, CA.
Designed for *everybody*.


@Therabody

Product Overview

SmartGoggles™

Relaxation meets comfort.

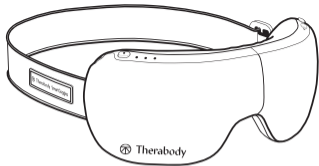From the creators of Theragun® comes SmartGoggles, a science-backed wearable that combines gentle heat, targeted vibration, and soothing massage with SmartSense Technology™, which uses your biometric data to personalize your experience and maximize your results. With three unique modes – SmartRelax, Focus and Sleep – SmartGoggles target pressure points to ease tension, relieve eye strain, ease discomfort related to headaches and migraines, and improve sleep quality. Pair SmartGoggles with TheraMind content in the Therabody App for an enhanced experience.

Intended Use

SmartGoggles intend to deliver comfort to the eyes and head to promote relaxation and sleep. Combining vibration with a gentle massage increases blood circulation to the eyes, reducing eye strain and irritation. Strengthen the connection between the mind and body with SmartGoggles.

EN Getting started

Step 1 - Download the Therabody app

**(If you've already used the Therabody app before, proceed to Step 2)*

1. Download the Therabody app from the App Store or Play Store.
2. Follow the in-app instructions to create a new app account and verify your email.
3. Once logged in, follow the onboarding steps to connect your SmartGoggles to the app. Be sure that the Bluetooth settings are enabled for your phone.
4. Select "SmartGoggles" as your Therabody device when prompted.
5. Long press the power button on your SmartGoggles to turn them on. SmartGoggles must be powered on to locate and pair to the app.
6. Tap "Pair Device" and select SmartGoggles to confirm the device and app connection.
7. Once connected, you can access preset routines to help you relax, focus, or sleep. You can customize the heat, vibration, and treatment length of these presets.
8. You can also select from a catalog of Theramind-enabled routines with science-backed audio to listen to during treatment for a fully immersed experience.

Step 2 - **(Start here if you've previously used the Therabody app)* Pair SmartGoggles with the Therabody app

- 4 1. Open the Therabody app and navigate to the Home screen.

2. Press the Bluetooth button in the lower right-hand corner. Tap “My Devices” to see a list of your Therabody devices. Tap the “+ Add New Device” button to connect your SmartGoggles. Ensure your phone’s Bluetooth settings are enabled.
3. Select “SmartGoggles” as your Therabody device when prompted.
4. Long press the power button on your SmartGoggles to turn them on. SmartGoggles must be powered on to locate and pair to the app.
5. Tap “Pair Device” and select SmartGoggles to confirm the device and app connection.
6. Once connected, you can access preset routines to help you relax, focus, or sleep. You can customize the heat, vibration, and treatment length of these presets.
7. You can also select from a catalog of Theramind-enabled routines with science-backed audio to listen to during treatment for a fully immersed experience.

Step 3 - Using your device

1. Place the device on clean, dry skin. Remove makeup, eyeglasses, and contact lenses prior to SmartGoggles use.
2. Place SmartGoggles on your head and adjust the strap. The Therabody SmartGoggles should perfectly align with your temples and sit comfortably on your nose bridge. The heart rate sensor should maintain contact with your cheekbone.
3. Keep your eyes closed during the treatment. Avoid frequent movement during treatment as it may impact heart rate tracking.
4. Long Press the Power Button to switch ON SmartGoggles.

5. Short Press the Power Button to toggle through the different modes, SmartRelax, Focus, and Sleep.
6. Short Press the Vibration Button to toggle through the four vibration patterns, low, high, wave, and pulse. Long Press the vibration button to stop the vibration.
7. Short Press the Heating Button to toggle through the three heating patterns. Long Press the heating button to stop the heating.
8. Long press the Power Button to turn off SmartGoggles. Make sure the Battery status LED is OFF.

Step 4 - Feel the difference.® Every day.

Whether used for alleviating eye strain or for general relaxation, SmartGoggles will help you feel your best.


See Warnings for important safety instructions.

Scan the QR code to get started.

Getting to know your wearable device

A SmartGoggles device

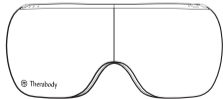B USB-C cable

C Soft case

D Instruction manual

Modes Overview

SmartGoggles has three strategically designed modes, SmartRelax, Focus, and Sleep.

Mode 1 - SmartRelax: SmartRelax promotes relaxation by gently massaging the eye area with heat. SmartRelax's vibration pattern reflects your heart rate, designed to optimize relaxation quickly.


Mode 2 - Focus: Focus aids in relaxation by releasing tension between the eyebrows. The single vibration pattern is intended to relax the mind and aid focus while heat is delivered to the eye area.

Mode 3 - Sleep: Sleep gently massages the temples and delivers slow-wave vibration to promote sleepiness.

SmartGoggles has two additional buttons for added customization of the three modes, heat and vibration, each with three intensity levels. See diagram on the next page for more details.

1. Four vibration settings (⚡): low, high, wave, pulse
2. Three heat settings (🔥): low, medium, high, off

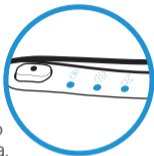P Power Button

Long press the power button to turn the device on and off.
Short press the power button while the device is powered off to check the battery status and to enable bluetooth.

Short press to toggle through Modes indicated by icons:

SmartRelax Mode

- V** **Vibration Settings**
Default: Pulse, Short press the vibration button to turn vibration off
- H** **Heat Settings**
Default: Low, Short press the heat button to toggle through additional heat settings: medium, high, off


Focus-Mode

- V** **Vibration Settings**
Default: Wave Short press the vibration button to toggle through additional vibration settings: low, high, off
- H** **Heat Settings**
Default: Medium, Short press the heat button to toggle through additional heat settings: high, off, low

Sleep Mode


- V** **Vibration Settings**
Default: Low Short press the vibration button to toggle through vibration additional settings: high, wave, off
- H** **Heat Settings**
Default: Off, Short press the heat button to toggle through additional heat settings: low, medium, high

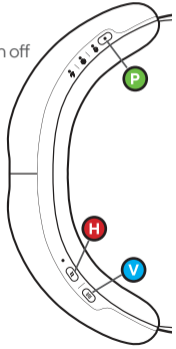Charging

- SmartGoggles are USB-C enabled
- SmartGoggles include a Therabody USB-C to USB-A cable.
- The USB port is located on the side of the device, behind the Power Button.
- Ensure that the battery status LED light is OFF before charging.
- Insert the USB-C cable into the port to charge SmartGoggles.
- The battery level is displayed in the LED indicator on the battery, next to the USB-C connector:
 - **Orange:** low charge, **Blue:** mid charge, **Green:** fully charged.
 - **When charging**, the LED will flash to indicate it is properly connected and will turn solid green when complete.

Note: If using a charger from a third party, ensure that it is from a certified manufacturer and has not suffered any structural damage. Do not try to use the device while charging.

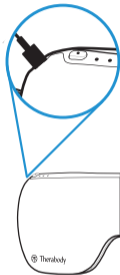Cleaning Instructions:

1. To clean, wipe the inner fabric of your SmartGoggles with a non-abrasive, lint-free cloth.
2. Lightly dampen a non-abrasive, lint-free cloth with water or an alcohol-based solution to clean if necessary.
3. Dry SmartGoggles with a non-abrasive, lint-free cloth.
4. SmartGoggles are NOT waterproof or machine washable. Do not allow excess liquids or water to enter the product.

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