ARTZT
vitality

WOBBLESMART®

MADE IN GERMANY

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ARTZT
vitality

The Wobblesmart® is used to improve body awareness and balance. Using the unstable platform for exercises requires stability and good coordination. It is mainly used for preventative conditioning, rehabilitation after injuries and in sports to increase performance.

Studies show that targeted body awareness exercises can reduce the frequency of falls and secondary injuries. By placing greater demands on postural muscles, the Wobblesmart® can also provide new energy to traditional exercises.

**How to use**

Make sure that the Wobblesmart® is resting firmly on the ground when exercising and that it cannot slip. Avoid training on slippery floors or rugs. If necessary, you can place the Wobblesmart® on an exercise mat.

Persons with significant disturbances of balance, posture or sensory function and persons with

epilepsy (due to medication, osteoporosis, low muscle strength) should only exercise under the guidance of a doctor, physiotherapist or specialised exercise instructor.


Remove sharp or pointed objects near the exercise area. Regular cleaning with soap and water is recommended. Use commercially available disinfectant spray for disinfection.

EXERCISE TIPS

Health-orientated, preventive exercises focus on improving general balance and body perception (proprioception). To achieve this goal, beginners should choose from the following exercises 2-4. The static holding exercises should be held 3 times for 10-30 seconds. The dynamic exercises should be repeated 10-15 times in 1-3 sets. At the very beginning of your exercises, you should pay particular attention to the slow and correct execution

of the exercises, even if they can only be held for 5 seconds at the beginning. Throughout the entire training, the strain should never feel excessive.

Once you have a certain foundation and you find the exercises easy, you can increase the number of exercises or exercise sequences, the difficulty level and the frequency.

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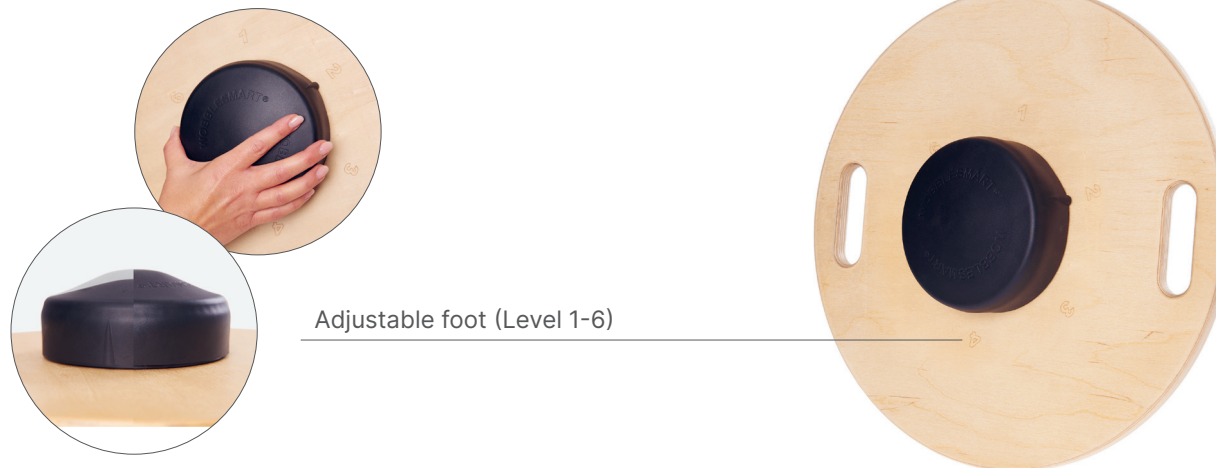Stand in the middle of the Woblesmart and find a stable position without wobbling (Woblesmart horizontal). Your arms should hang loosely beside your body or be held slightly apart.

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Slowly turn your head to the right and across the centre to the left. If this goes well, try it faster without losing your stable stance.

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Move your head up and down (chin touches chest). Vary the range of movement and speed without losing your stable stance. Exercises A-C can also be done with your eyes closed to increase the difficulty.

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Stand in the middle of the Wobblesmart and find a stable position without wobbling (Wobblesmart horizontal). Your arms should hang loosely beside your body.

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Slowly and purposefully shift your weight forwards and backwards (Wobblesmart tilts forwards and backwards) and to the side (Wobblesmart tilts to the right and left),...

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Stand with your feet very close together in the centre of the Wobblesmart, stretch your arms out to the sides and find a stable position.

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Raise your bent left leg forward and return to the starting position. Change legs.

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Stretch your left leg out to the side and return to the starting position. Change legs.

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Stand in the middle of the Wobblesmart and find a stable position without wobbling Wobblesmart horizontal). Your arms should hang loosely beside your body.

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Bend your knees and bring your stretched arms forward parallel to your bent knees. Slowly return to the starting position.

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Like B but with a loop above the knees. Choose a wide position for your feet so the loop is under tension. Maintain this tension throughout the exercise.

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
Stand behind the Wobblesmart. Place your right foot in the centre of the Wobblesmart and shift your weight to the front leg. Slowly release your left foot from the floor and find a stable standing position on your right leg.

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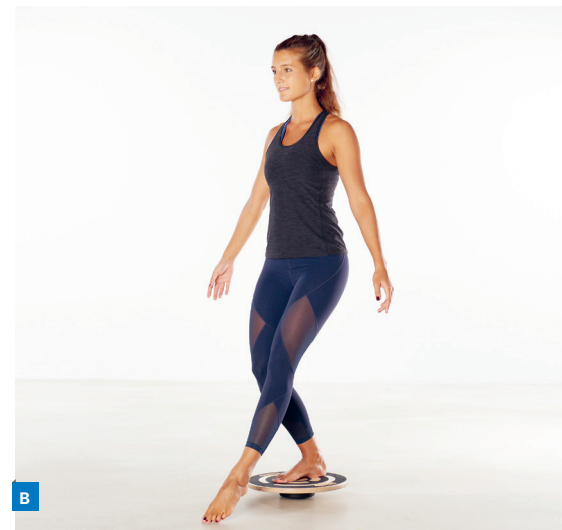With the standing leg slightly bent, first touch the floor at '12 o'clock' with the big toe of the left leg and then touch the lower leg.

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Then move from 11 to 6 o'clock. Change legs. (Move your left leg from 12 to 6 o'clock and do the same with the right leg.)

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Stand behind the Wobblesmart. Place your right foot in the centre of the Wobblesmart.

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Shift your weight onto your right foot and carefully lift your left foot off the floor. Find a stable position. Change legs.

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If you do the preparatory exercise well, you can do a dynamic 'running movement' with your moving leg and arms in a one-legged stance. Try to maintain a stable standing position throughout the movement.

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Place both palms shoulder-width apart and flat on the Woblesmart and get into a press-up position. Find and hold a stable support position with the Woblesmart aligned parallel to the floor.

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If exercise A is successful, grasp the recessed grips of the Woblesmart and return to the press-up position with your legs extended.

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Perform press-ups and try to keep the Woblesmart parallel and your elbows as close to your torso as possible. Look towards the ground.

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Place both palms shoulder-width apart and flat on the Wobblesmart and get into a press-up position. Find a stable support position with the Wobblesmart aligned parallel to the floor.

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Now alternately lift one leg off the floor and keep your pelvis parallel to the floor.

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Get into a press-up position with the Wobblesmart under the balls of your feet. Place your hands under your shoulders and have your feet open and hip-width apart. Keep your pelvis parallel to the floor. You can also do this exercise with forearm support.

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When you can hold the position well, slowly rotate your pelvis from left to right (shifting your weight from left to right foot). Try to keep your shoulder girdle stable.

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Lie on your back with your knees bent and the Wobblesmart under the soles of your feet. Your arms should rest alongside your body.

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Raise and lower your pelvis in a controlled and slow manner.

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Sit in the centre of the Wobblersmart and lift both feet off the floor. Your knees should be bent. Your hands should rest loosely on your opposite shoulders. Find a stable position.

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Open your arms quickly and rhythmically, close them again and find a stable sitting position as quickly as possible.

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Stand behind the Wobblesmart. In a wide stride position, place your right foot on the centre of the Wobblesmart. Shift your weight centrally between your legs and rest on the ball of the rear foot.

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Now bend and straighten both knees and keep the Wobblesmart as parallel to the floor as possible and your upper body upright. Change legs.

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Attach the Super Band textile to the left hand grip and stand on the Wobblesmart with your feet far apart. Grasp the band loop with both hands (fist) in front of your left hip.

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Now pull the band slowly and in a controlled manner diagonally upwards to the right. Switch sides.

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All rights reserved. Printed in Germany.
© 2022 Ludwig Artzt GmbH, Schiesheck 5,
Ludwig Artzt GmbH - Schiesheck 5 - 65599 Dornburg - Germany