DEUTSCH

Wichtige Information

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ENGLISH

Important Notice

Dear Customer, Congratulations on your purchase of a high-quality AIREX brand product. To enjoy the benefits of our CE-conforming AIREX products as long as possible, please observe these care instructions.

Application The AIREX Balance products are ideal for coordination and balance training. Prior to use, the training surface must be checked to ensure it is slip-resistant. A non-slip mat (e.g. AIREX mat) must be used if applicable. Wearing shoes can reduce the lifespan of AIREX Balance products.

Appearance Due to the special manufacturing process, certain deviations or imperfections such as small surface bubbles, minor size differences or slight discolouration may occur. This does not influence the functionality in any way.

Characteristics AIREX soft foams do not absorb water and feature permanent buoyancy in water. Pressure loss in the individual cells occurs when pressure loads are applied for long periods of time. While this results in pressure marks and creases, the material usually reforms gradually. Special hygienic properties offer permanent, effective protection against bacteria and fungi.

Cleaning The preferred cleaning method is to use a soft brush and lukewarm water with detergent. If disinfection is required, a conventional disinfectant such as Merfen, Descogen or chromium may be used. Treatment with benzene or spot remover as well as hot water or steam sterilisation (autoclave sterilisation) is not permissible.

Storage AIREX Balance products are best stored lying flat. To avoid pressure marks, do not leave objects resting on the material for extended periods of time. Continuous pressure loads for more than approximately 12 hours should be avoided. Bending in storage leaves visible pressure marks and creases where the material does not fully reform. Ongoing storage in water and extended exposure to direct sunlight reduce the lifespan of the products, since the surface may become brittle and the material stiff over time. To extend the lifespan, we recommend rinsing the product with clear water after use, shaking off the water and allowing the product to air dry. For hygienic reasons, damp products should not be stored in lockers or stacked.

We hope you enjoy and have great success with the products you have purchased. Please contact your specialist dealer if you have further questions.

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Manufactured by:
Airex AG
Industrie Nord
5643 Sins, Switzerland
www.my-airex.com


<p>Produktions-Los-Nummer: Production lot number:</p>
<p>Always indicate for further enquiries or complaints.</p>

www.my-airex.com | facebook.com/myairex | made in switzerland

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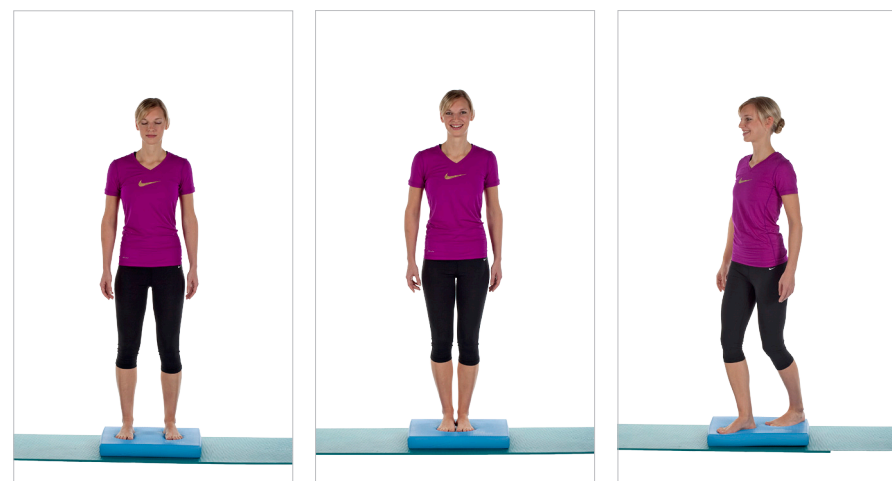
1. Standing on two legs
Stand with feet hip-width apart in erect position. Control your balance on your feet without bending your hip or helping with your arms.
Variation: Close your eyes, vary position of the feet.

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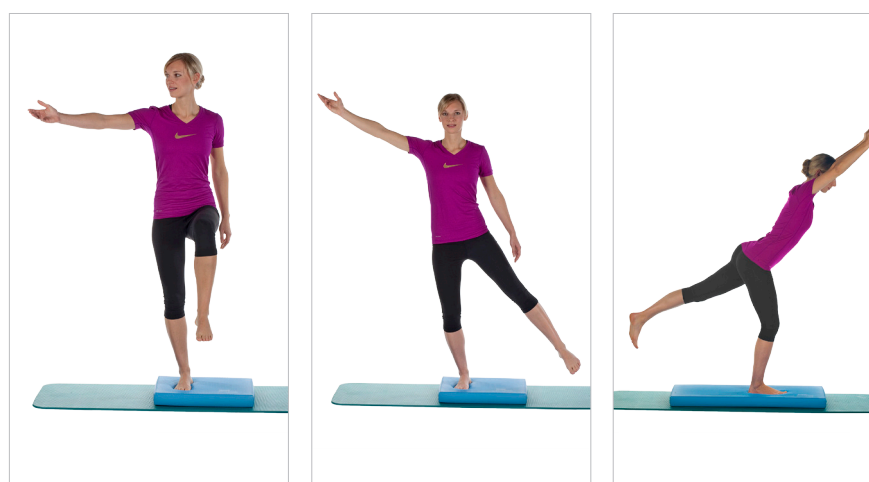
2. Standing on one leg
Stand on one leg, flex the free leg at the hip and lift straight up, extend the opposite arm away from the body, look straight ahead.
Variation: Change arms and legs.

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3. Running imitation
Stand on one leg, imitate a running movement with your arm and leg.
Variation: Change the dynamics of bending and stretching.

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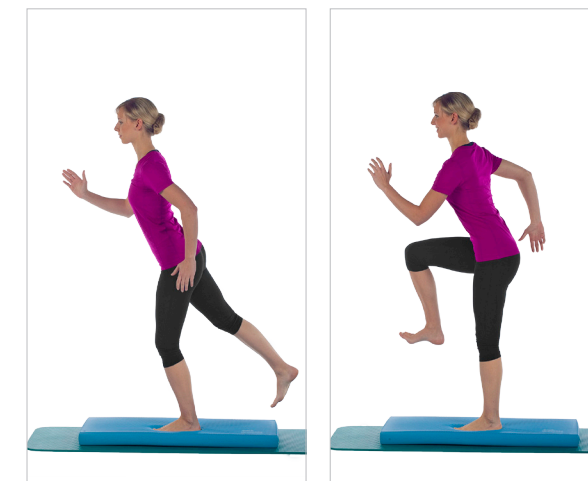Warm-up: Before starting to exercise, a 10-minute warm-up is recommended – e.g. walking / running on the spot or on the AIREX Balance-pad Elite.

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Exercises: Take care to carry out the movements in a slow and controlled way. Start with a few repeats and slowly increase their number. Later several series are possible.

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Quality beats quantity! Diversify your workout! Alternate your exercises regularly. Work out as often as possible. All exercises can be performed with the Balance-pad Elite as well as the Balance-pad XLarge. Use an AIREX mat as base for the exercises.

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Cool Down: Enjoy a few moments of relaxation. Lie on your back, legs bent, place hands on your stomach and close your eyes. Breathe in and out deeply and evenly.

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8. Thoracic spine extension
Lie on your back on the Balance-pad, arms in extension of the body.
Variation: Also extend your legs.

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

7. Plank exercise
Support yourself on your forearms on the Balance-pad and stretch your hip. Stabilize your core.
Variation: Lift knees from the floor, prop up on your hands.

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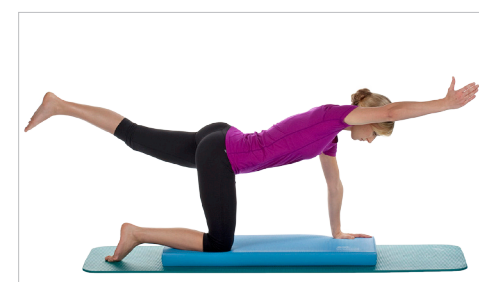
6. Tabletop
Assume tabletop position and touch your knee with the opposite elbow.
Variation: Change arms and / or close your eyes.

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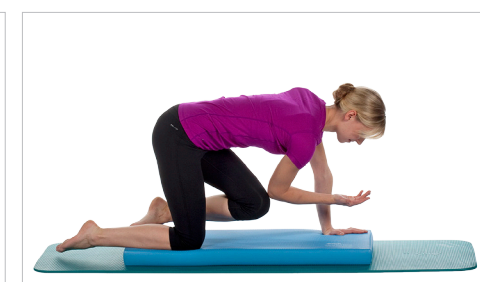
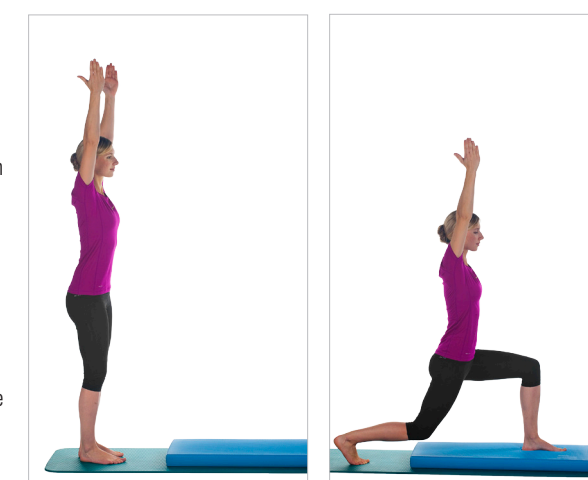
4. Forward lunge
Stand in an upright position in front of the Balance-pad and lunge forwards. Return to the original position, ensuring your core is upright and your arms are stretched.
Variation: Change legs.

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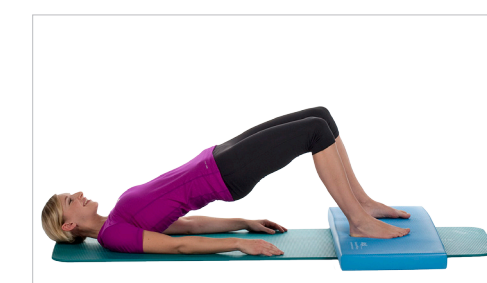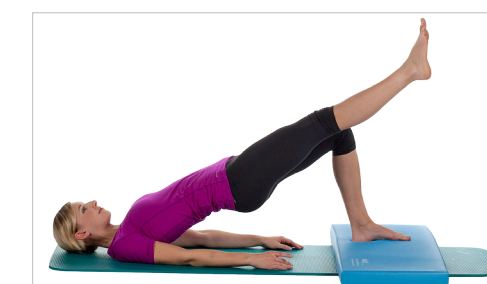5. Pelvis lift
Lie on your back and slowly raise and lower your pelvis.
Variation: Alternately lift your feet with a stretched leg and lower them back onto the Balance-pad.

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